


Stone Retrieval Basket

The Stone Retrieval basket is designed to manipulate, grasp and remove calculi and foreign bodies from the urinary tract. Its detachable handle allows the withdrawal and re-entrance of the scope when the basket is situated.

■ Features

- Torque-resistant Nitinol wire basket with round construction facilitating efficient capture and atraumatic manipulation of calculi
- Low friction kink-resistant sheath providing smooth endoscopic passage
- Available with different designs: tipless designs allowing for closer proximity to stones; Tipped baskets facilitating access to stones in torturous anatomy
- Ergonomic handle promoting excellent single-hand control through the procedure
- Handles and sheaths can be detached easily and completely when needed; the nitinol core can be left in place working as a guide wire when the basket cannot be released from the stone; the basket can be easily re-assembled once it is removed from the patient; Not suitable for 3-wire stone retrieval basket

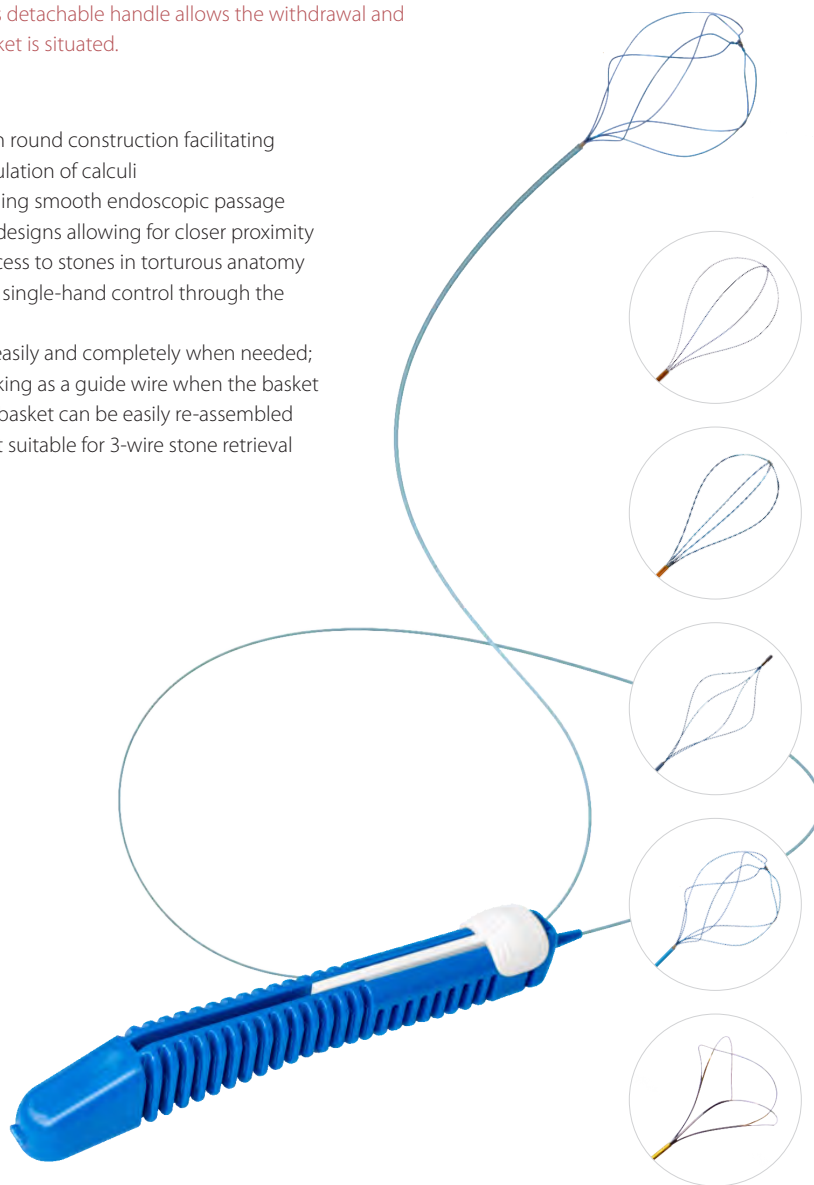

Size (Fr)	Working Length (cm)	Basket Description	Basket OD (mm)	Cat. No.
1.8	90	Straight 4 wire, zero tip	11	U04B010110
1.8	120	Straight 4 wire, zero tip	11	U04B010210
2.2	120	Straight 4 wire, zero tip	10	U04B010310
2.5	120	Helical 4 wire, straight tip	15	U04B010410
2.5	120	Helical 6 wire, zero tip	17	U04B010510
1.7	130	3 wire, open tip	8	U04B011010
1.7	130	3 wire, open tip	11	U04B011110

* 1pcs/box, 10pcs/ctn